

Stream, Continuous Filter Holder



WATCH VIDEO

121°C / 1bar / 20 min



Stream 3 / Stream 10



Stream, Continuous Filtration System

Features

- ▶ A smart kit enables continuous sample feeding. No more repetitive pouring tasks.
- ▶ Equipped with a sinker that misses no drops till the end.
- ▶ Filtration adaptor can be connected to all GL45 bottles.

Applications


Large-volume liquid purification and filtration, such as solvents used in HPLC, GC, AA analysis, and culture medium purification.

Specification

Material	Stream 3	Stream 10
Continuous filtration kit	PTFE, Silicone	PTFE, Silicone
Funnel	Borosilicate glass	Borosilicate glass
Support base	Borosilicate glass	Borosilicate glass
Membrane support	Sintered glass	Sintered glass
Clamp	Aluminum	Aluminum
Filtration adaptor	PP	PP

Specification

Funnel capacity	300 mL	1000 mL
Filter diameter	47 mm / 50 mm	90 mm
Effective filtration area	13.2 cm ²	45.5 cm ²
Hose barb	ID4 ~ ID8	ID4 ~ ID8

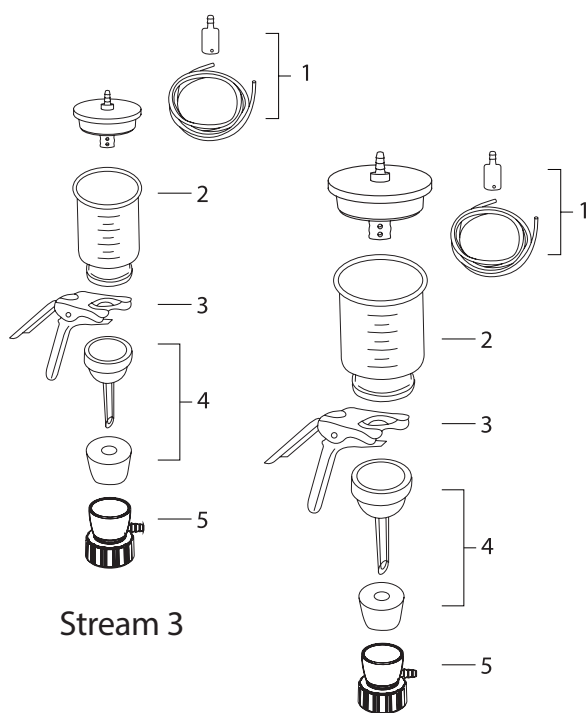

Stream 3

Stream 10

Ordering Information

Stream 3, Continuous Filter Holder (300 mL)

167200-S3

1	Continuous filtration kit for 300 mL glass filter holder	167210-36
2	Glass funnel, 300 mL	167210-03
3	Aluminum clamp, 47 mm	167240-01
4	Sintered glass support base, 47 mm (Including No.8 silicone stopper)	167230-03
5	GL45 filtration adaptor	197000-65

Stream 10, Continuous Filter Holder (1000 mL)

167200-S10

1	Continuous filtration kit for 1000 mL glass filter holder	167210-35
2	Glass funnel, 1000 mL	167210-10
3	Aluminum clamp, 90 mm	167240-02
4	Sintered glass support base, 90 mm (Including No.8 silicone stopper)	167230-11
5	GL45 filtration adaptor	197000-65