

Large Number of Samples

Alligator X MultiVac Vacuum Filtration System

Efficient / Space-Saving / Versatile



Compliance with ISO 8199


Features



Save Time

1

- ◆ Built-in liquid pump transfers filtered liquids directly to waste.
- ◆ Quick-Fit hose barb can be quickly connected to either side of the manifold towards the vacuum source.

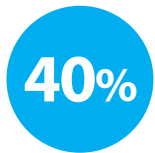

Save 40% Bench Space

2

- ◆ No need of waste bottle saving much bench space.



Versatile

3

- ◆ Replaceable funnel adaptor makes it compatible with various kinds of filter holders.



Certification

CE / RoHS / IP 64




▲ Compatible with various kinds of filter holders.

Applications

- ◆ Microbiological test by Membrane Filter Technique
- ◆ Suspended solid test for water monitoring

Specification

Model	Alligator 200	
Power	DC24V, 1.8A	
Max. vacuum	313 mbar	
Max. liquid flow rate	4.5 L/min	
Motor	60W Brushless DC Motor, improves the pump life and performance	
Noise level	60 dB	
Hose barb	ID8	
Material contact with filtrate	SS316, FKM, Silicone, PP, PPS	
International certification	CE, RoHS, IP 64	
Dimension (LxWxH)	23.5 x 19 x 20.2 cm	
Model	MultiVac 301 / 601 - MB	MultiVac 310 / 610 - MS
Number of branch	3 / 6	3 / 6
Filter diameter	47 mm	Optional
Hose barb	ID8	ID8
Material	SS316	Aluminum

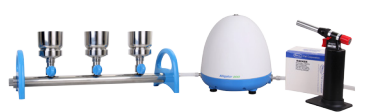

MultiVac 301 - MB - A, Vacuum Filtration System



MultiVac 310 - MS - A, Vacuum Filtration System



MultiVac 601 - MB - A, Vacuum Filtration System



MultiVac 610 - MS - A, Vacuum Filtration System



181200-01(02)	Alligator 200, Diaphragm Liquid Pump with AC100~240V Adaptor, US Plug (EU Plug)
180301-01(02)-A	MultiVac 301-MB-A, Vacuum Filtration System with AC100~240V Adaptor, US Plug (EU Plug)
180601-01(02)-A	MultiVac 601-MB-A, Vacuum Filtration System with AC100~240V Adaptor, US Plug (EU Plug)
180311-01(02)-A	MultiVac 310-MS-A, Vacuum Filtration System with AC100~240V Adaptor, US Plug (EU Plug)
180611-01(02)-A	MultiVac 610-MS-A, Vacuum Filtration System with AC100~240V Adaptor, US Plug (EU Plug)
180100-01	SF Stainless Steel Funnel 100 ml
180100-03	SF Stainless Steel Funnel 300 ml
180100-05	SF Stainless Steel Funnel 500 ml
200300-01	MF 3, Magnetic Filter Holder 300 ml
200500-01	MF 5, Magnetic Filter Holder 500 ml
190100-73	SS316 Adaptor for Millipore disposable plastic funnel
190100-74	SS316 Adaptor for Sartorius disposable plastic funnel
190100-75	SS316 Adaptor for Pall disposable plastic funnel

Rocker Scientific Co., Ltd.

Headquarter

19F.-6, No.206, Guanghua 1st Rd., Lingya Dist.,
Kaohsiung City 802, Taiwan (R.O.C.)
T/ +886-7-227-3358 F/ +886-7-227-3716

Taipei office

11F., No.6, Zhulin Rd., Linkou Dist.,
New Taipei City 244, Taiwan (R.O.C.)
T/ +886-2-2603-3311 F/ +886-2-2603-6622

Get more information for Rocker

www.rocker.com.tw
export@rocker.com.tw

