

BioDolphin Suction Kit



BioDolphin Suction Kit

Introduction

BioDolphin is a versatile aspiration tool that can convert any vacuum source to aspirate liquids from virtually any laboratory container, such as cell culture dishes, 96-well plates, flasks, centrifuge tubes, bottles, and so on. The ergonomic handle with one touch and lock switch is simply designed for fitting with different aspiration adapters. 1-/8-channel tip adaptor with ejector are designed to simplify routine work, and unique 15cm tip adaptor is dually suitable for 200µl and 1000µl tips. The whole BioDolphin suction kit except storage rack is autoclavable to prevent cross contamination and keep your safety from liquid biological waste.

Benefits

- **Versatile**
 Various adapters allow to aspirate liquids from all kinds of containers
- **Labour-Saving**
 Single-handed tip ejection fittings simplify tip discarding procedure
- **Ergonomic**
 Handle and lock switch are convenient and effective for long-time or repetitive working
- **Easy-to-use**
 Easy plug design helps quick change for all fittings
- **Autoclavable**
 All fittings except rack are autoclavable which can prevent any contamination

Applications


- Cell culture aspiration
- Laboratory waste suction
- ELISA wash solution suction

Video Introduction



Standard accessory

1. Handle
2. 1-channel tip adaptor (50mm)
3. 1-channel SS adaptor (80mm)
4. 1-channel SS adaptor (40mm)
5. 8-channel SS adaptor (40mm)
6. 8-channel tip adaptor with ejector
7. Rack
8. 1-channel tip adaptor with ejector
9. 1-channel tip adaptor (150mm)

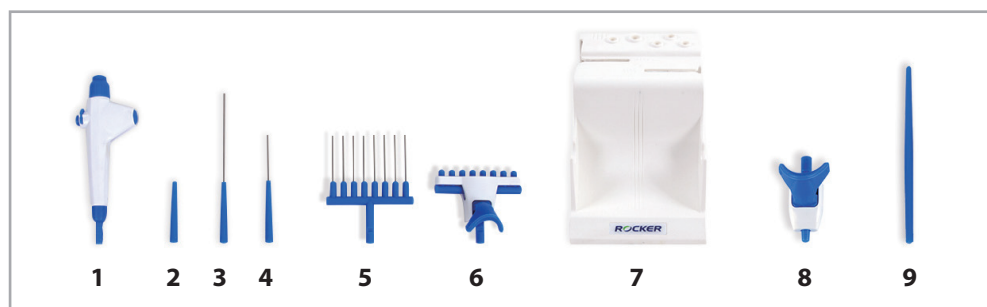

Product portfolio

According to laboratory requirements, BioDolphin can be used with various types of vacuum pump or suction system.



Lafil 100



Lafil 200 Plus



Lafil 200-BioDolphin



199100-00	BioDolphin - LB Suction Kit (Light Blue)
197100-01	Lafil 100, Portable Suction System with AC100-240V adaptor, US plug
197100-02	Lafil 100, Portable Suction System with AC100-240V adaptor, EU plug
197201-01	Lafil 200 Plus, Portable Suction System with AC100V-240V adaptor, US plug
197201-02	Lafil 200 Plus, Portable Suction System with AC100V-240V adaptor, EU plug
197202-01	Lafil 200-BioDolphin, Portable Suction System with AC100V-240V adaptor, US plug
197202-02	Lafil 200-BioDolphin, Portable Suction System with AC100V-240V adaptor, EU plug
199100-30	Handle
199100-46	1-channel tip adaptor (50mm)
199100-47	1-channel S.S adaptor (40mm)
199100-48	1-channel S.S adaptor (80mm)
199100-49	8-channel S.S adaptor (40mm)
199100-50	8-channel tip adaptor with ejector
199100-61-1	Rack
199100-75	1-channel tip adaptor with ejector
199100-76	1-channel tip adaptor (150mm)

ROCKER

11F., No.6, Zhulin Rd., Linkou Dist.,
New Taipei City 244, Taiwan (R.O.C.)

T/ +886-2-26033311

F/ +886-2-26036622

<http://www.rocker.com.tw>

Distributor