

Filtration Accessories



GL45 Filtration Adaptor

Specification	GL45 filtration adaptor
Catalog number	197000-65
Material of cap	PP
Material of gasket	PTFE
Stopper	No.8 stopper
Hose barb	ID4 ~ ID8
Autoclavable	Yes



GL45 Suction Cap NEW

Specification	GL45 suction cap	SC45 Suction Cap	SC45 Suction Cap Kit
Catalog number	197000-60	197000-61-P2N	197000-61-P2
Material of cap	PP	PP	PP
Material of gasket / O-ring	Silicone	EPDM	EPDM
Overflow protection	-	-	Yes
Hose barb	ID6	ID4 ~ ID8	ID4 ~ ID8
Autoclavable	Yes	Yes	Yes

SC45 Suction Cap GL45 Bottle Applications



► On waste bottle

Collecting liquid waste by being connected to vacuum pumps and filtration accessories.



► On safety bottle

Working as a safety bottle to prevent liquid from being sucked into pump.




► On moisture trap bottle

Attaching to GL45 bottle filled with desiccant to prevent moisture from entering the vacuum pump.




► Vapor absorption (neutralization) bottle

Ideal for absorbing acid or alkaline vapor to prevent harmful vapors from spreading into the environment or pump.

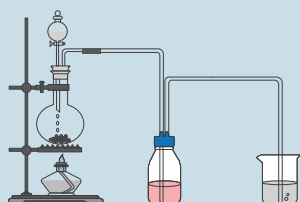

► Transfer bottle

Ideal for transferring sterile media stored in GL 45 bottle to another container or reactor with peristaltic pump in a closed system.

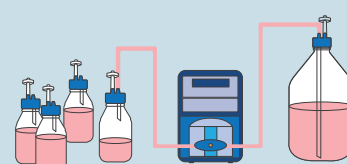

* Please refer to P.92 to select laboratory glass bottle.